RT-MIS

ENGLISH	2
FRANCAIS.....	23
DEUTSCH.....	44
ESPAÑOL	66
ITALIANO	87
PORTUGUES	108

ENGLISH

GENEXPATH RT-MIS software instruction for use

User precautions.

***In vitro* diagnostic medical device according to Directive (EU) 98/79/EC**

For *in vitro* diagnostic use

It is for professional use only.

Read all information in this user guide before use.

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Telephone : +41 (044) 805 76 81

Fax : +41 (044) 805 76 75

This document and its content are the property of GENEXPATH. They are only intended for contractual use by customers for the product described herein, excluding any other use. This document and its content may not be used or distributed for any other purpose, nor shared, disclosed or reproduced in any way without the prior written permission of GENEXPATH. GENEXPATH does not grant any license through this document.

The instructions given in this document must be strictly and explicitly followed by qualified and trained personnel to ensure appropriate use of the product described herein. All content in this document must be thoroughly read and understood before using the described product.

© 2025 Genexpath. All rights reserved.

ENGLISH	2
Foreword - please read carefully	5
About this guide	5
Product description	5
General information	5
GENEXPATH LymphoSign	5
GENEXPATH SarcomaFusion	5
GENEXPATH CustomPanel	6
System requirements	6
Home page	7
Analyses	9
Start a new analysis	9
Results	10
Processing details	11
All tests	11
LymphoSign Test	13
SarcomaFusion test	17
Specificities	19
Performances	19
Repeatability	19
LymphoSign Test	19
Random forest	19
Analytical specificity	19
Analytical sensitivity	19
Graph representation	19
SarcomaFusion Test	20
Filtering non-specific fusion	20
Display of the fusions	20
Associated Bibliography	21
Internal control	21
Correction of hopping index	21
Notes	22
Reference	22

Foreword - please read carefully

GENEXPATH RT-MIS is a comprehensive bioinformatics solution which includes different data processing algorithms. It should be used critically by professionals in human genetics. Although GENEXPATH strives to offer high-quality software, it cannot guarantee the accuracy of the information and predictions it provides.

About this guide

This guide describes how to use the software GENEXPATH RT-MIS, Version 06, published in October 2025.

Product description

General information

After completing the *in vitro* protocol of a Genexpath test and sequencing, the software **GENEXPATH RT-MIS** allows the loading of sequencer files in FASTQ format, as well as the list of barcodes used during testing.

This software carries out demultiplexing to assign sequences to each sample. It then precisely identifies the gene expression markers and quantifies them.

It assesses the sequencing quality of each sample by quantifying the number of reads identified and the number of UMI (unique molecular identifiers) detected.

For each sample, it generates a report to help the user analyse the data. This data can be downloaded.

GENEXPATH LymphoSign

Through the use of LymphoSign tests, the software helps analyse results obtained after completing the *in vitro* protocol and sequencing. The analysis is based both on quantification of gene marker expression, but also the presence or absence of mutations and chromosome translocations.

The artificial intelligence classification model was trained with a database containing over 3000 analyses.

If this database were to change, users will be kept informed via the RT-MIS platform.

The expression level of each marker is normalised in relation to the average gene expression of all markers studied, allowing a comparison of analyses.

GENEXPATH SarcomaFusion

Through the use of SarcomaFusion tests, the software helps analyse results obtained after completing the *in vitro* protocol and sequencing. RT-MIS provides information on the presence or absence of chromosome translocations amongst 210 determined fusions associated with

sarcomas. It also establishes a link with bibliography references associated with detected fusions.

The bibliography used was created from Mitelman databases [1] and PubMed research for all fusion transcripts associated with sarcomas. It is non-exhaustive and specific to sarcomas.

GENEXPATH CustomPanel

In the context of the use of CustomPanel assays, the software allows the analysis of the results obtained after the *in vitro* protocol and sequencing have been performed. The raw data and expression level of each marker are provided. The expression level of each marker is normalized in relation to the average gene expression of all the markers studied, thus allowing the comparison of the analyses with each other.

System requirements

The software RT-MIS is directly accessible through the user space at <https://connect.genexpath.com>. Access is restricted to user accounts authorised by Genexpath.

To get an account, send an email to the following address: support@genexpath.com.

An internet connection is required. The software is compatible with the browsers Mozilla Firefox, Microsoft Edge and Google Chrome. It does not require the installation of additional tools.

Home page

After logging into the user area, access to RT-MIS is provided by clicking the **RT-MIS** button in the left-hand menu. A dedicated page for the platform will open, containing several sections:

- **Presentation:** briefly summarizes how RT-MIS works and includes a link to the support email.
- **Online Documentation:** provides links to Genexpath product inserts.
- **Latest News:** displays the most recent information about updates.

By clicking the "**ACCESS MY RT-MIS ACCOUNT**" button (Figure 1), you will be redirected to your RT-MIS user account.

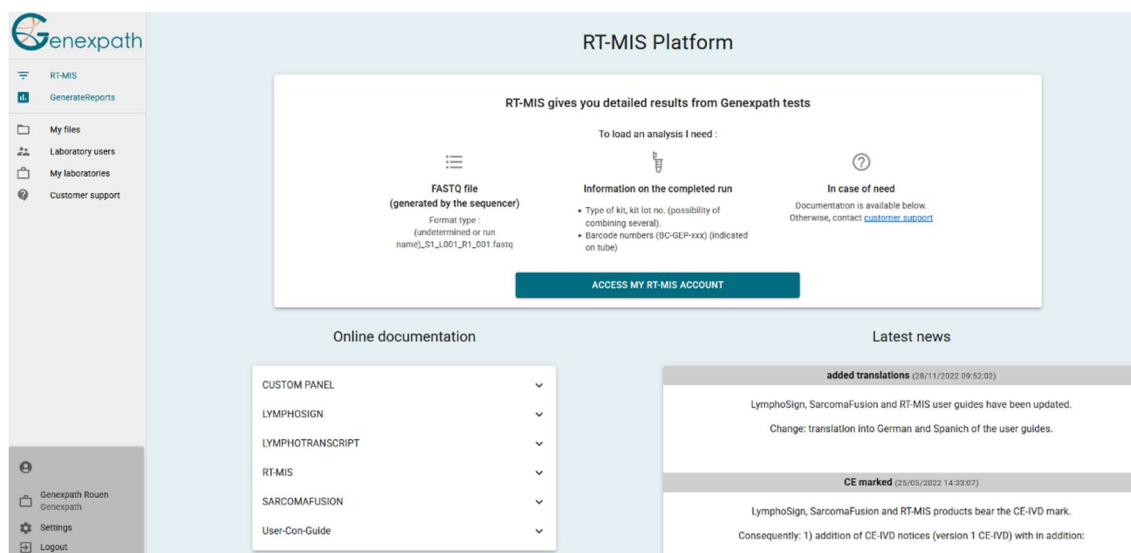
Figure 1: User space page

The home page is shown (Figure 2):

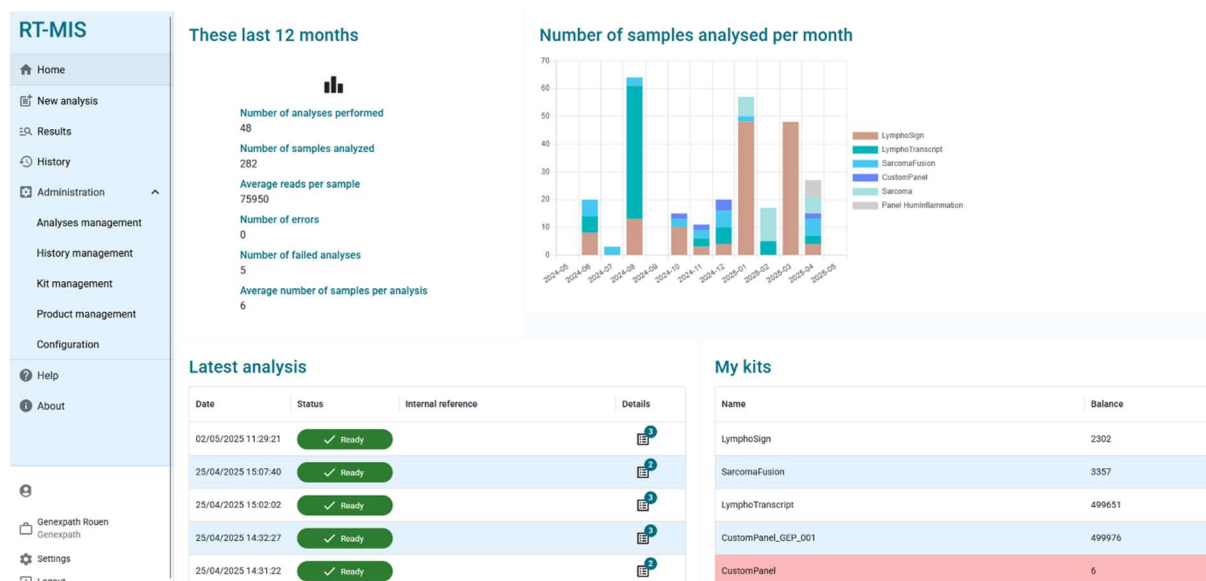
Figure 2: RT-MIS home page

It contains several sections:

- **"Number of Samples Analyzed Per Month"**: a graph showing the number of samples analyzed, grouped by kit type.
- **"These Last 12 Months"**: a summary of statistics (including the number of analyses performed, samples analyzed, errors, failed analyses, and the average reads per sample and samples per analysis).
- **"Latest Analyses"**: displays the most recent analyses performed, their status, and access to detailed results by clicking the corresponding icon. The number shown in the top right corner indicates the number of samples processed (and therefore deducted from the corresponding kit balance).
- **"My Kits"**: shows information about the kits assigned to the laboratory the user is affiliated with, including the number of remaining analyses. If the balance falls below 7, the background turns red.

On the left, the main menu (top) allows you to access:

- On the home page
- New analysis launch
- Display of results
- kits 'history
- Interactive help (for the home page and the "new analysis" page)
- About


The main menu (bottom) indicates the login with which you are connected and allows:

- To return to the user area homepage or the RT-MIS homepage
- To choose to receive an email at the end of an analysis
- Sign out.

Analyses

Start a new analysis

- 1) In the menu on the left, click on **“New analysis”** (Figure 3).

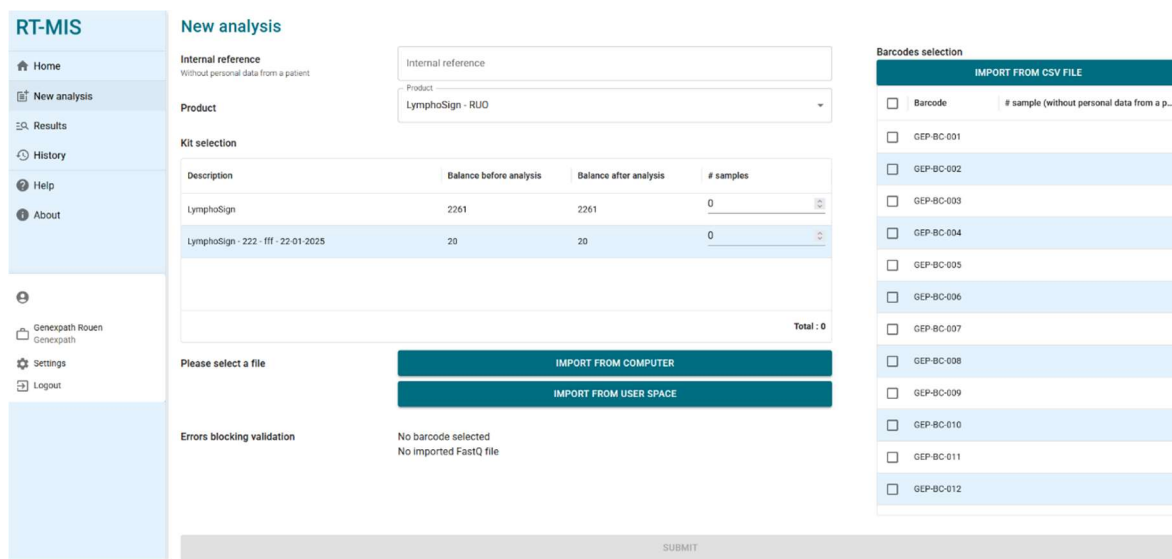RT-MIS

New analysis

Internal reference
Without personal data from a patient

Product
LymphoSign - RUO

Kit selection

Description	Balance before analysis	Balance after analysis	# samples
LymphoSign	2261	2261	0
LymphoSign - 222 - If - 22-01-2025	20	20	0
			Total : 0

Please select a file

IMPORT FROM COMPUTER

IMPORT FROM USER SPACE

Errors blocking validation

No barcode selected
No imported FastQ file

Barcodes selection

IMPORT FROM CSV FILE

☐ Barcode # sample (without personal data from a p...

☐ GEP-BC-001

☐ GEP-BC-002

☐ GEP-BC-003

☐ GEP-BC-004

☐ GEP-BC-005

☐ GEP-BC-006

☐ GEP-BC-007

☐ GEP-BC-008

☐ GEP-BC-009

☐ GEP-BC-010

☐ GEP-BC-011

☐ GEP-BC-012

SUBMIT

Figure 3: “New analysis” page

- 2) Indicate an internal reference to the analysis with any patient data information (figure 4)
- 3) Select in the dropdown menu “product”, the used product.
- 4) In “kit selection”, indicate the number of analyses that you need to count down per kit (according to the batch number used in the “*in vitro*” part.
- 5) Select a FastQ file to analyse by clicking on:
 - **“IMPORT FROM COMPUTER”** for a file saved on the computer, or
 - **“IMPORT FROM USER SPACE”** for a file already imported to the user space.

The FASTQ file must respect the naming convention of Illumina sequencers.

MiSeq example: X_S1_L001_R1_001.fastq.gz

The software checks the format and integrity of the FastQ file(s) before proceeding to the analysis. This check lasts from a few seconds to a few minutes, depending on the file size.

- 6) Tick (right side of the screen) the barcodes used for the analysis (numbers printed on the tubes).

For each selected barcode, an analysis will be deducted.

You can add an identifier for each selected barcode.


No patient data can be entered.

You can also have the option of downloading directly a CSV file with its information from the CSV file example.

Issues that could block the validation are:

- Absence of FASTQ file
- Barcodes selection mismatch

7) Click on **“SUBMIT”**.

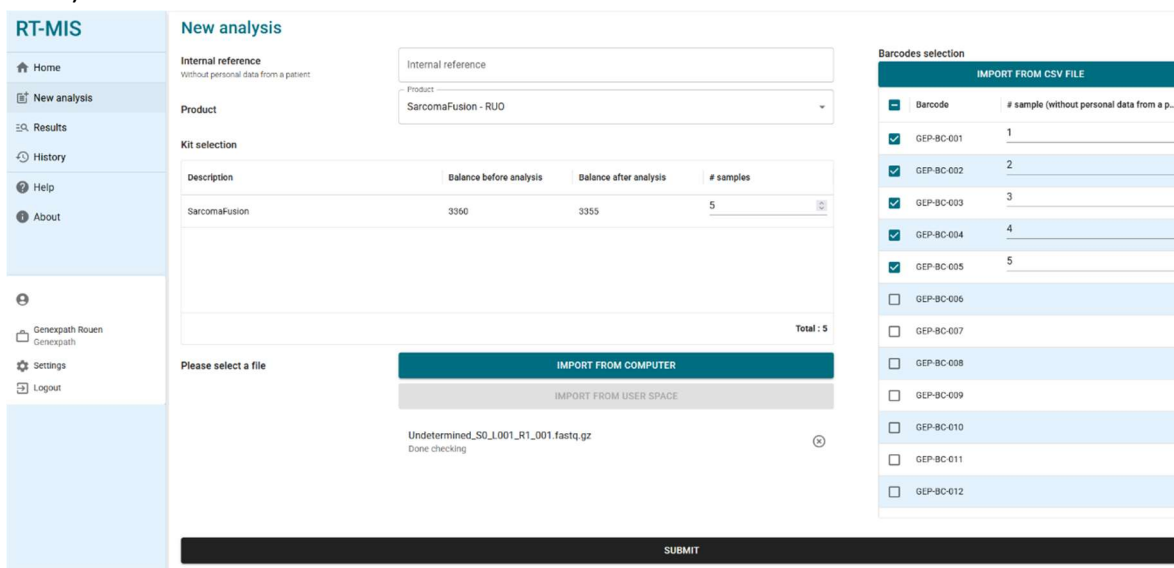RT-MIS

New analysis

Internal reference
Without personal data from a patient

Internal reference

Product
SarcomaFusion - RUO

Kit selection

Description	Balance before analysis	Balance after analysis	# samples
SarcomaFusion	3360	3355	5

Total : 5

Please select a file

IMPORT FROM COMPUTER

IMPORT FROM USER SPACE

Undetermined_S0_L001_R1_001.fastq.gz
Done checking

Barcodes selection

IMPORT FROM CSV FILE

Barcode	# sample (without personal data from a p...
☒ GEP-BC-001	1
☒ GEP-BC-002	2
☒ GEP-BC-003	3
☒ GEP-BC-004	4
☒ GEP-BC-005	5
☐ GEP-BC-006	
☐ GEP-BC-007	
☐ GEP-BC-008	
☐ GEP-BC-009	
☐ GEP-BC-010	
☐ GEP-BC-011	
☐ GEP-BC-012	

SUBMIT

Figure 4: Information on the user analysis

RT-MIS opens the **“Results”** page.

Results

On the left side menu, click on item **“Results”** (figure 5).

This page offers access to new and previous results.

RT-MIS indicates the analysis status: **“pending processing”**, **“processing”** or **“processing complete”**.

The results will be available once processing is complete. Processing can last several minutes.

To access analysis details, click on the corresponding icon .

RT-MIS

Home

New analysis

Results

History

Administration

- Analyses management
- History management
- Kit management
- Product management
- Configuration

Help

About

Julien Binet

jbinet@genexpath.com

Genexpath Rouen

Genexpath

Settings

Results

Date	Product	Internal reference	User	Status	Details
02/05/2025 11:29:21	LymphoSign		Laurine BECQUET	✓ Ready	3
25/04/2025 15:07:40	CustomPanel		Laurine BECQUET	✓ Ready	2
25/04/2025 15:02:02	LymphoTranscript - RUO		Laurine BECQUET	✓ Ready	3
25/04/2025 14:32:27	SarcomaFusi		Julien Binet	✓ Ready	3
25/04/2025 14:31:22	LymphoSign		Julien Binet	✓ Ready	3
08/04/2025 10:55:16	Panel Huminflar		Laurine BECQUET	✓ Ready	1
02/04/2025 11:35:19	SarcomaFusion		Julien Binet	✓ Ready	3
02/04/2025 11:34:12	LymphoSign		Julien Binet	✓ Ready	3
01/04/2025 14:41:42	Sarcoma V2		Laurine BECQUET	✓ Ready	3
17/03/2025 09:54:09	LymphoSign		Laurine BECQUET	✓ Ready	44
04/03/2025 14:38:50	LymphoSign		Laurine BECQUET	✓ Ready	1
17/02/2025 14:35:10	LymphoTranscript - RUO		Laurine BECQUET	✓ Ready	3
12/02/2025 10:18:44	Sarcoma		Laurine BECQUET	✓ Ready	3
05/02/2025 13:48:56	Sarcoma		Laurine BECQUET	✓ Ready	3

Total Rows: 218

Figure 5: "Results" page

Processing details

All tests

This page offers a summary of the analysis characteristics (Figure 6).


- The status
- The centre
- The project
- The internal reference
- The product used
- The user
- The FASTQ file save date
- The FASTQ file processing date
- The kit used and the number of analyses used on it.

Process detail

Status	✓ Ready	
Center	Genexpath	
Laboratory	Genexpath Rouen	
Internal reference		
Product	LymphoSign - RUO	
User		
Record date	05/08/2024 13:19:39	
Process date	05/08/2024 13:20:02	
Kit	Kit	Analysis consumed on the kit
	LymphoSign	11

Figure 6: "Processing details" page

If any of the barcodes that are populated are not found in the parsed FASTQ file, a warning message will be displayed. You can request a re-load directly by clicking on the corresponding button.

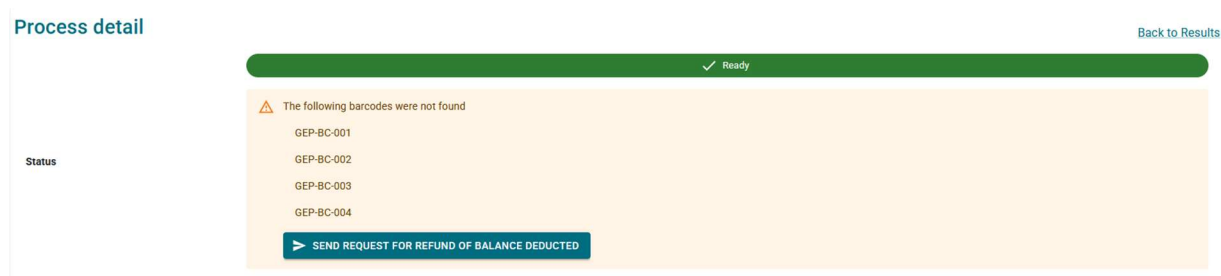Quality metrics are represented for the entire treatment (Figure 7):

- Genexpath barcodes: Percentage of reads associated with Genexpath libraries on the entire FASTQ file loaded.
- Percentages of barcoded reads according to its barcode
- Representation of the reads of the loaded FASTQ file

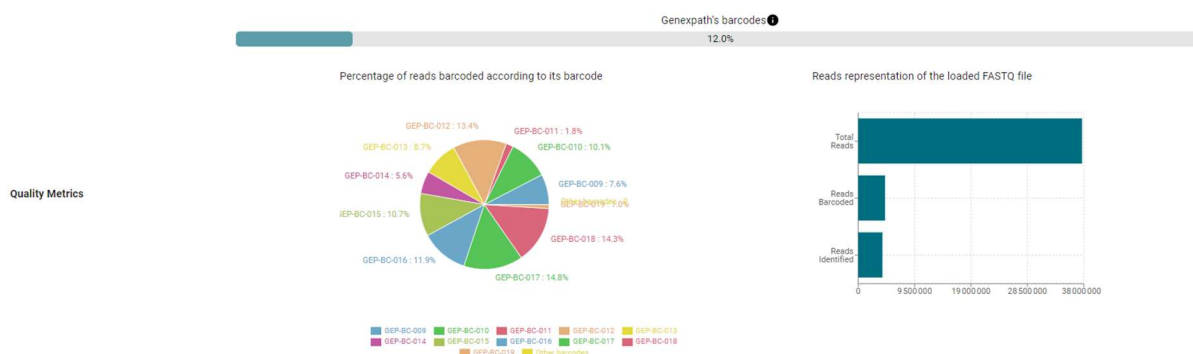
Figure 7: Quality metrics

For each sample tested, you will find the results of your analysis with a reminder of the barcode used and the sample reference, if applicable. This part allows you to observe the informativity for all the samples in the analysis (Figure 8).

	Barcode	Sample	Informativity LymphoSign	Reports	Count data
Barcodes	 CEP-BC-014		91.9%		
	 CEP-BC-015		90%		

Figure 8: Access to analyses by barcode

LymphoSign Test

It is possible to generate a customizable analysis report by clicking on  .

A report generation page opens, by default on tab “results”:

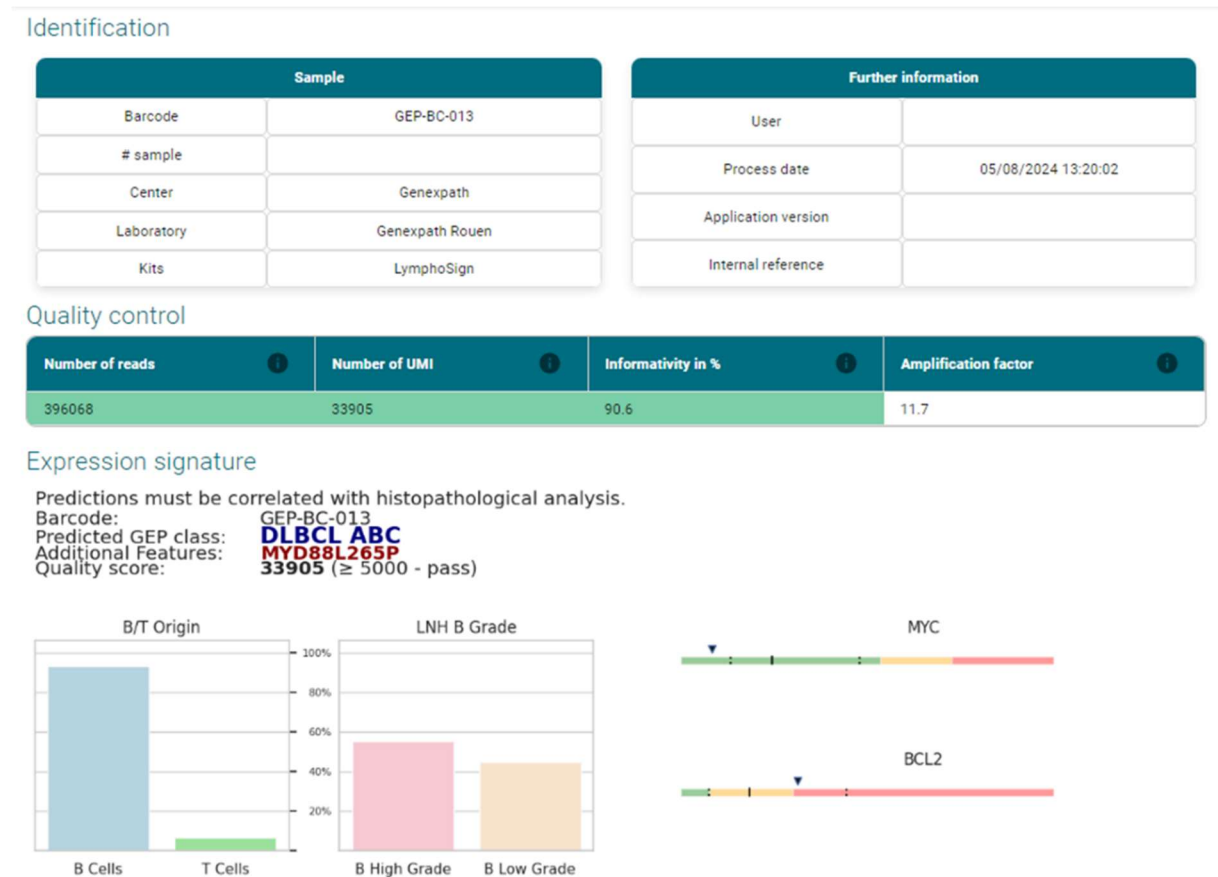
Figure 9: Example report for LymphoSign test

By default, the report contains the following information:

- Analysis identification: reminder of analysis characteristics
- Quality control: number of reads and UMI (Unique Molecular Identifiers), the informativity (number of junction L probes – R waited probes / total number of reads *100) and the amplification factor.
- The expression signature: barcode, predicted class, optional additional feature, quality score, as well as a set of graphs detailing the inferences made to predict the class of non-Hodgkin's lymphoma.

Name	Description
Barcode	Reminder of barcode used for the analysis
Predicted GEP class	Class based on AI-predicted gene expression profile
Additional Features	Additional information
Quality score	Quality score based on number of UMI

You can also choose which items you will eventually see appear in the report using the "Add Item" buttons depending on how many charts you will be able to add. In addition, you can view the report before downloading by clicking on "view PDF". The graph is downloaded by clicking on "download PDF".

The additional elements are:

- A customizable graph (Figure 10) representing the level of expression of several genetic markers of interest according to the position of your sample (red dot) as well as:
 - o the entire database
 - o One of the two classes LNH_B or LNH_T
 - o Of a sub-entity

Custom graph

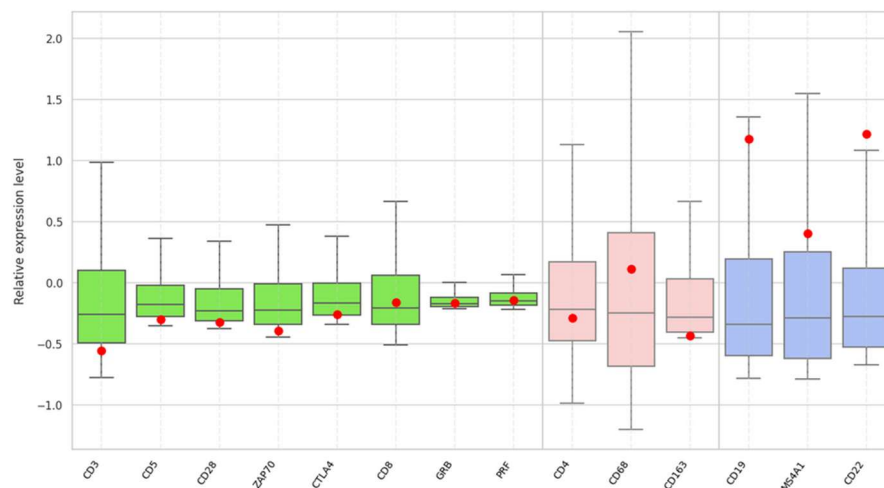
Figure 10: Generation of custom graph

You can at the top right (Figure 11):

- o Choose the number of classes wished by clicking on "add a class"
- o Determine the name of each class
- o Choose the markers of interest by clicking on the pencil of the corresponding class

Edit class

Class name

CD3

CD5

CD28

ZAP70


CTLA4

CD8

GRB

PRF

Associated color

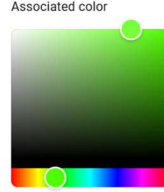SUBMIT

Figure 11: Determination of markers for a given class

- Predefined graphical representations with the comparison of the level of expression (Figure 12). The reference expression levels are displayed according to the chosen class and the sample analyzed is represented by a black diamond or a red dot depending on the tab. It is possible to compare gene expression levels.

Tab	Description
Comparison	Comparison of expression profiles between two LymphoSign signature classes (the predicted one and another chosen by the user)
B vs T	Specific B and T lymphocyte markers
IG	Different immunoglobulins
MUTS	Specific markers such as point mutations
DLBCL	Specific markers of diffuse large B cell lymphomas (DLBCL)
PTCL	Specific markers of peripheral T-cell lymphomas
Other	Specific markers of different lymphomas
Prognosis	Markers which could affect patient prognosis

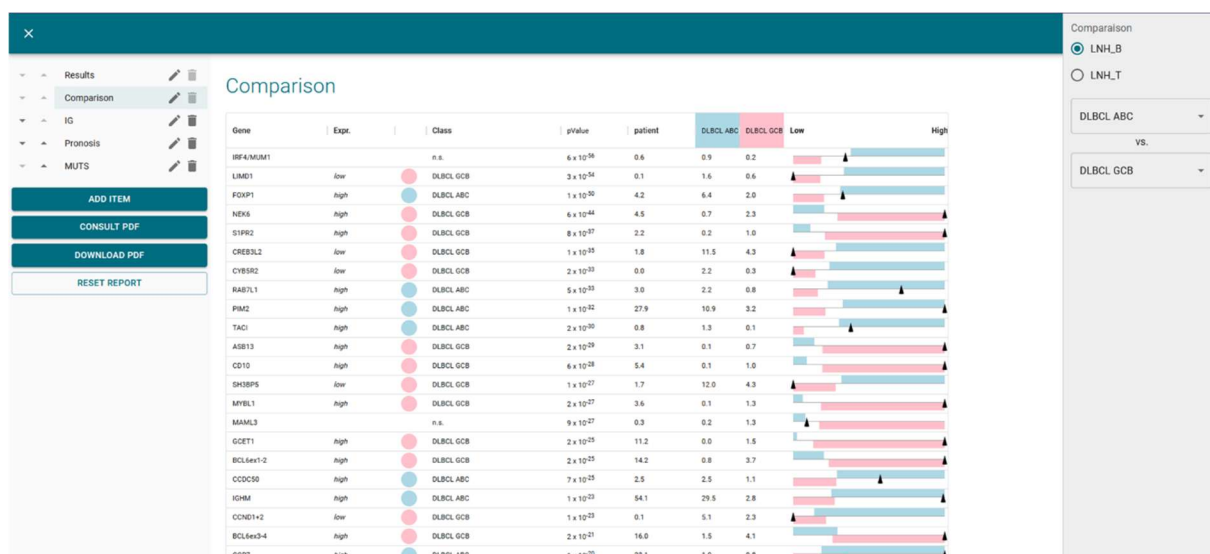
Figure 12: Example of advanced interpretation with comparison between PMBL and GCB

Only the first 30 genes of the two-class comparison are displayed.

Each of these comparisons can be personalised using the right-hand menu, through which you can select the type of classes to compare (Figure 13):

- The entire database
- one of the two classes LNH_B or LNH_T
- A sub-entity

Only the first 30 genes of the comparison of 2 classes will be display.

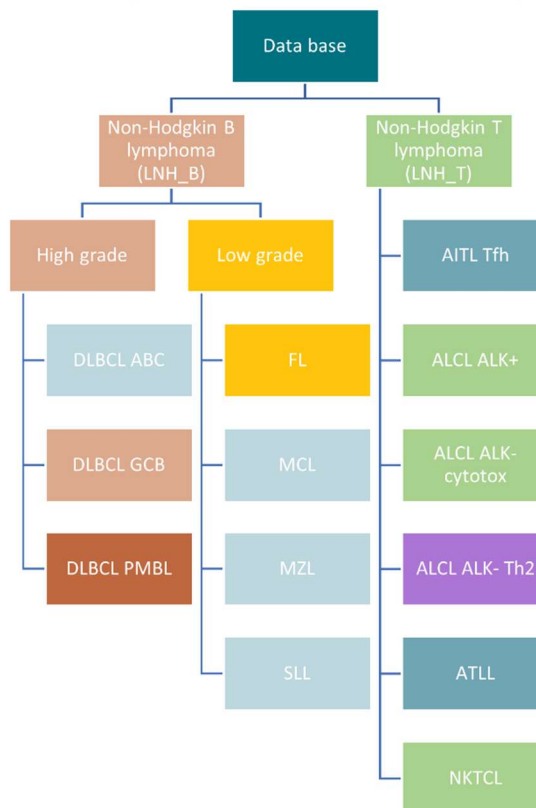
Figure 13: Representation of different classes predicted by RT-MIS

You may also download in “**management details**” page:

- For each sample: the CVS file of all raw and relative data of by clicking on 
- The CSV file with all raw and relative data by clicking on “**DOWNLOAD ALL COUNTING DATA**”
- All pdf report or CSV files can be downloaded by clicking on “**ALL REPORT**”

SarcomaFusion test

By clicking on  icon of barcode of interest, you can also generate a report (figure 14) with or without the associated bibliography, consult and download it.

Identification

Sample		Further information	
Barcode	GEP-BC-002	User	
# sample		Process date	01/08/2024 13:37:01
Center	Genexpath	Application version	
Laboratory	Genexpath Rouen	Internal reference	
Kits	SarcomaFusion		

Quality control

Quality control SarcomaFusion		Quality control GAPDH	
Number of reads	189638	Number of reads	187109
Number of UMI	2954	Number of UMI	11767
Amplification factor	64.2	Amplification factor	15.9
Informativity		Informativity GAPDH	
85.4%		95.7%	

Results

Fusion	Number of reads	Number of UMI	Amplification factor	References
✓ PAX3_7E7G-FOXO2D	158329	1163	136.1	7 references

Fusion : PAX3_7E7G-FOXO2D

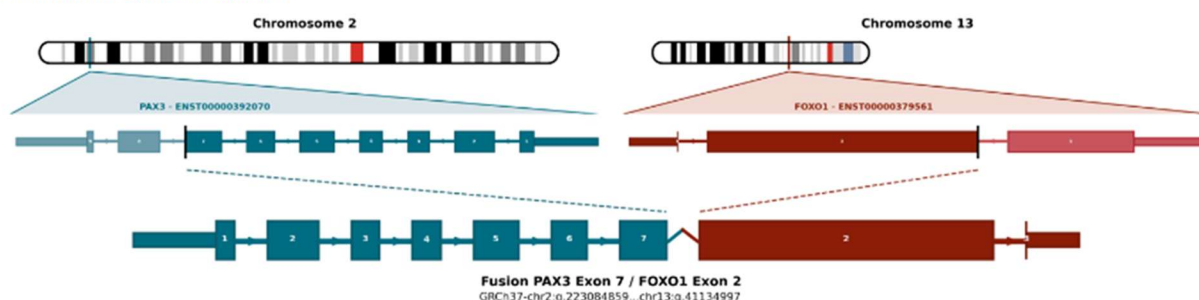
Figure 14: "Report" page for the SarcomaFusion test

Information available:

- Result with
 - Summary of the analysis characteristics
 - Quality control
 - number of reads and UMI for the sample of interest.
 - The informativity in % as well as the amplification factor.
 - number of reads and UMI for the GAPDH internal control
 - The informativity in % as well as the amplification factor for internal control

- The results of the fusion may have been recognized
- The bibliography references if the fusion is associated to a sarcoma
- The raw data obtained by RT-MIS, without filtration of the number of UMI, with a graphic representation of the detected fusion depending on the number of reads or UMI and a summarize table.
- A graphic representation of the raw data

You may also download in “**management details**” page:

- For each sample: the CVS file of all raw and relative data of by clicking on 
- The CSV file with all raw and relative data by clicking on “**DOWNLOAD ALL COUNT DATA**”
- All pdf report or CSV files can be downloaded by clicking on “**DOWNLOAD ALL REPORT**”

Fusion transcript notation in RT-MIS is as follows (Figure 15):

Nomenclature : JAZF1E3G-SUZ12E2D

Figure 15: Schematic representation of naming system for fusion transcripts

By clicking on “n references”, the bibliography details open in a new page (Figure 16).

Bibliographic report

Name	References
JAZF1/SUZ12	Novel EPC1 gene fusions in endometrial stromal sarcoma. Genes Chromosomes Cancer
JAZF1/SUZ12	Cytogenetic and molecular genetic analyses of endometrial stromal sarcoma: nonrandom involvement of chromosome arms 6p and 7p and confirmation of JAZF1/JJAZ1 gene fusion in t(7;17). Cancer Genet Cytogenet
JAZF1/SUZ12	Molecular analysis of the JAZF1-JJAZ1 gene fusion by RT-PCR and fluorescence in situ hybridization in endometrial stromal neoplasms. Am J Surg Pathol
JAZF1/SUZ12	Extrauterine endometrial stromal sarcoma with JAZF1/JJAZ1 fusion confirmed by RT-PCR and interphase FISH presenting as an inguinal tumor. Virchows Arch
JAZF1/SUZ12	A neoplastic gene fusion mimics trans-splicing of RNAs in normal human cells. Science
JAZF1/SUZ12	Frequent fusion of the JAZF1 and JJAZ1 genes in endometrial stromal tumors. Proc Natl Acad Sci U S A

Figure 16: “Bibliography references” page

Specificities

Performances

Repeatability

The repeatability of the RT-MIS software is defined as its ability to obtain the same results. The same FASTQ analyzed in triplicate with 15 samples for each of two tests was studied. The RT-MIS software has a repeatability of 100%. The same FASTQ file containing several samples gives:

- for a LymphoSign test: the same classification as well as the same expression levels for each of the markers
- for a SarcomaFusion test: the same result (presence or absence of a fusion) as well as the associated bibliography.

LymphoSign Test

Random forest

Normalization of expression data is a step aimed at making samples sequenced under different experimental conditions comparable (quantity and quality of RNA, type of sequencer, number of reads generated, etc.). In order to erase these experimental biases, the raw counting data is divided by the median of the counts of a fixed list of markers.

From the normalized expression data, RT-MIS classifies the sample by Random Forest (RF). This classification is based on the application to the expression data of 10 RF classification models from the Machine Learning software from the database of nearly 3000 cases. The different probabilities of these classifiers make it possible to identify different signatures and classes described in Figure 13.

Analytical specificity

A study was carried out on a larger scale with 283 cases (Bobee et al., Blood Cancer Journal, 2020). The objective of this study was to demonstrate the relevance/power of the Random Forest algorithm allowing classification by comparing the results with another classification technique. 94.3% of the results were concordant. Thanks to the algorithm, some of the cases were later reclassified.

Analytical sensitivity

The artificial intelligence was trained on a database of nearly 3,000 cases. The prediction is optimal when there are at least 5000 UMI.

Graph representation

The level of expression of each of the markers corresponds to the raw count data divided by the median of the counts of a fixed list of markers.

The data is represented in the form of boxplots. They correspond for a given marker to the reference values of the database according to the entire database, of one of the two classes LNH_B or LNH_T or of a sub-entity.

Each boxplot is organized as follows:

- The central value of the graph is the median
- The edges of the rectangle are the quartiles (Q1 at the bottom and Q3 at the top)
- The endpoints of the boxplot are calculated using 1.5 times the interquartile space (the distance between the 1st and 3rd quartiles).

SarcomaFusion Test

Filtering non-specific fusion

The background noise is established by a list of non-specific fusions which is filtered to facilitate the understanding of the results. It was created from the results obtained from 5 tested samples whose expected fusions were known. Non-specific fusions were detected in at least 2 samples.

This list is as follows:

SS18E10G-SS18E3D	FUSE5G-STAT6E17D	RAD51BE3G-TFE3E7D	PHF1E12G-DDIT3E3D
NAB2E3G-DDIT3E3D	HEY1E4G-CREB3L2E6D	FUSE9bG-CAMTA1E9D	PHF1E12G-BCORE8D
PAX3_7E7G-TFE3E3D	KIF5BE24G-CAMTA1E9D	CTNNB1E1G-NTRK1E16D	PHF1E12G-RETE8D
ZC3H7BE12G-RELA1E11D	RBPMSE5G-EGFE17D	COL1A1E37G-STAT6E17D	FUSE5G-ALKE13D
TEAD1E8G-ETV1E4D	EWSR1E15G-ALKE20delD	CTNNB1E1G-NTRK3E18D	FUSE5G-ALKE16D
CITED2E2G-ERGE10D	LIFRE1G-NFATC2E3D	COL1A1E37G-STAT6E17D	COL1A1E47G-TFE3E3D
HEY1E4G-NTRK1E10D	TFGE4G-NTRK3E11D	CTNNB1E1G-NTRK3E18D	VCLE16G-KLF17E2D
TFGE7G-FLI1E8D			

These non-specific fusions do not appear in the PDF report, nor in the CSV file which can be downloaded.

Display of the fusions

The pairs of probes detected correspond to potential genetic anomalies. They are compared to the bibliographic database. The objective is to determine whether the fusion detected has already been described in the literature in the context of sarcomas or not. Depending on this referencing, the fusion is retained or filtered according to the algorithm below (Figure 17):

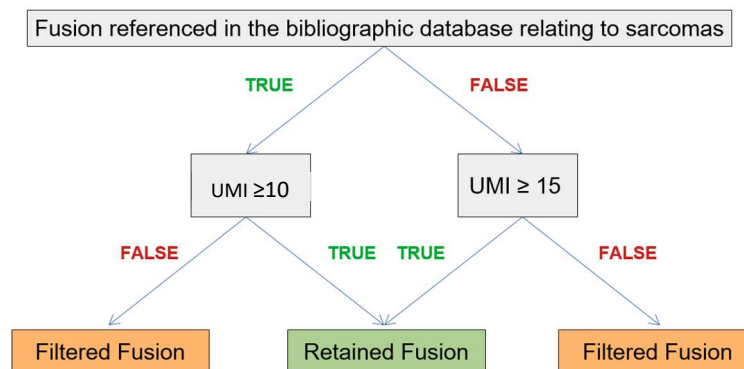
Figure 17: Algorithm of filtration representation

The retained fusion (or fusions) is returned as a result in the "Filtered data" part of the report and, if applicable, associated with its literature and the filtered fusions are entered in the "Raw data" part.

Associated Bibliography

The bibliographic database used was created from the Mitelman databases [1] and Pubmed searches for all the fusion transcripts associated with sarcomas. It is non-exhaustive and specific to sarcomas. It will be updated twice a year.

Internal control

It is recommended to perform an internal control during the *in vitro* test using the same barcode and control barcode number for the same sample, as indicated in the SarcomaFusion instructions for use. During computer analysis, RT-MIS combines the raw data from both barcodes to generate a single report.

The sample with the barcode must be sequenced as a normal sample with 100 000 reads. Nevertheless, its usage has no impact of the analysis countdown.

Correction of hopping index

Hopping index is a phenomenon, during sequencing, of exchanges of barcodes leading to a wrong attribution of the sequence/sample couple. This phenomenon, in the case of the SarcomaFusion tests, can lead to misinterpretation of the results.

To correct this technical effect, the UMIs assigned to each pair of left/right probes of the different sequenced samples are compared with each other. The UMI sequences shared by the same pair of probes between two barcodes of the same run are eliminated when they are not amplified.

Notes

The company **GENEXPATH** does not permanently store results generated by the software **GENEXPATH RT-MIS**. The data must be downloaded directly after each analysis and saved by the user in their document management system.

However, analysis results are available for one year.

Information in this document is likely to change.

Reference

[1] Mitelman Database of Chromosome Aberrations and Gene Fusions in Cancer (2021). Mitelman F, Johansson B and Mertens F (Eds.), <https://mitelmandatabase.isb-cgc.org>

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RT-MIS	
Nom	Référence
PAX3/FOXO1	Alveolar rhabdomyosarcoma of the head and neck region in older adults: genetic characterization and a review of the literature. Hum Pathol
PAX3/FOXO1	Cryptic insertion of 3FOXO1 into inverted chromosome arm 2q in the presence of two normal chromosome 13s and 13 small interstitial duplications in a patient with alveolar rhabdomyosarcoma. Cancer Genet
PAX3/FOXO1	Molecular Cytogenetics Detect an Unbalanced t(2;13)(q36;q14) and PAX3-FOXO1 Fusion in Rhabdomyosarcoma With Mixed Embryonal/Alveolar Features. Pediatr Blood Cancer
PAX3/FOXO1	Molecular diagnostics in the management of rhabdomyosarcoma. Expert Rev Mol Diagn
PAX3/FOXO1	Clinicopathologic and Molecular Features of a Series of 41 Biphenotypic Sinonasal Sarcomas Expanding Their Molecular Spectrum. Am J Surg Pathol
PAX3/FOXO1	Early detection of the PAX3-FOXO1 fusion gene in circulating tumor-derived DNA in a case of alveolar rhabdomyosarcoma. Genes Chromosomes Cancer
PAX3/FOXO1	Alternate PAX3-FOXO1 oncogenic fusion in biphenotypic sinonasal sarcoma. Genes Chromosomes Cancer

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Process detail

[Back to Results](#)

Status	✓ Ready	
	⚠ The following barcodes were not found GEP-BC-001 GEP-BC-002 GEP-BC-003 GEP-BC-004 ➤ SEND REQUEST FOR REFUND OF BALANCE DEDUCTED	

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Identification

Sample		Further information	
Barcode	GEP-BC-013	User	
# sample		Process date	05/08/2024 13:20:02
Center	Genexpath	Application version	
Laboratory	Genexpath Rouen	Internal reference	
Kits	LymphoSign		

Quality control

Number of reads	Number of UMI	Informativity in %	Amplification factor
396068	33905	90.6	11.7

Expression signature

Predictions must be correlated with histopathological analysis.

Barcode: GEP-BC-013
 Predicted GEP class: **DLBCL ABC**
 Additional Features: **MYD88L265P**
 Quality score: **33905** (≥ 5000 - pass)

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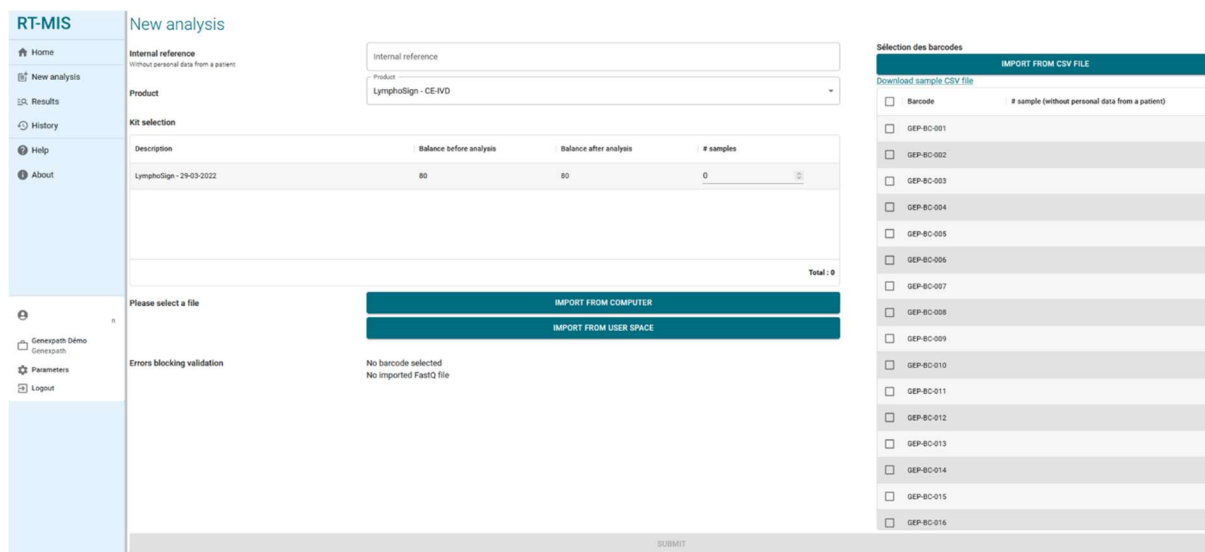
The screenshot displays the 'New analysis' page. On the left is a sidebar with navigation links: Home, New analysis (selected), Results, History, Help, and About. Below these are links for Genexpath Demo, Genexpath, Parameters, and Logout. The main content area is titled 'New analysis' and includes a 'Please select a file' section with two buttons: 'IMPORT FROM COMPUTER' and 'IMPORT FROM USER SPACE'. Below these buttons, it states 'Errors blocking validation' and 'No barcode selected. No imported FastQ file'. On the right, there is a 'Selection des barcodes' panel with a table of barcodes (GEP-BC-001 to GEP-BC-016) and checkboxes for selection. A 'SUBMIT' button is at the bottom right.

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RT-MIS		Results				
Home		Date	Product	Internal reference	User	Status
New analysis		02/05/2025 11:29:21	LymphoSign		Laurine BECQUET	✓ Ready
Results		25/04/2025 15:07:40	CustomPanel		Laurine BECQUET	✓ Ready
History		25/04/2025 15:02:02	LymphoTranscript - RUO		Laurine BECQUET	✓ Ready
Administration		25/04/2025 14:32:27	SarcomaFusi		Julien Binet	✓ Ready
Analyses management		25/04/2025 14:31:22	LymphoSign		Julien Binet	✓ Ready
History management		08/04/2025 10:55:16	Panel Huminfan		Laurine BECQUET	✓ Ready
Kit management		02/04/2025 11:35:19	SarcomaFusion		Julien Binet	✓ Ready
Product management		02/04/2025 11:34:12	LymphoSign		Julien Binet	✓ Ready
Configuration		01/04/2025 14:41:42	Sarcoma V2		Laurine BECQUET	✓ Ready
Help		17/03/2025 09:54:09	LymphoSign		Laurine BECQUET	✓ Ready
About		04/03/2025 14:38:50	LymphoSign		Laurine BECQUET	✓ Ready
		17/02/2025 14:35:10	LymphoTranscript - RUO		Laurine BECQUET	✓ Ready
		12/02/2025 10:18:44	Sarcoma		Laurine BECQUET	✓ Ready
		05/02/2025 13:48:56	Sarcoma		Laurine BECQUET	✓ Ready

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Process detail					
Status	✓ Ready				
Center	Genexpath				
Laboratory	Genexpath Rouen				
Internal reference					
Product	LymphoSign - RUO				
User					
Record date	05/08/2024 13:19:39				
Process date	05/08/2024 13:20:02				
Kit	<table> <tr> <th>Kit</th><th>Analysis consumed on the kit</th></tr> <tr> <td>LymphoSign</td><td>11</td></tr> </table>	Kit	Analysis consumed on the kit	LymphoSign	11
Kit	Analysis consumed on the kit				
LymphoSign	11				

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Back to Results

✓ Ready

△ The following barcodes were not found

GEP-BC-001

GEP-BC-002

GEP-BC-003

GEP-BC-004

➤ SEND REQUEST FOR REFUND OF BALANCE DEDUCTED

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RT-MIS		Results				
Home		Date	Product	Internal reference	User	Status
New analysis		02/05/2025 11:29:21	LymphoSign		Laurine BECQUET	✓ Ready
Results		25/04/2025 15:07:40	CustomPanel		Laurine BECQUET	✓ Ready
History		25/04/2025 15:02:02	LymphoTranscript - RUO		Laurine BECQUET	✓ Ready
Administration		25/04/2025 14:32:27	SarcomaFusi		Julien Binet	✓ Ready
Analyses management		25/04/2025 14:31:22	LymphoSign		Julien Binet	✓ Ready
History management		08/04/2025 10:55:16	Panel HumInfar		Laurine BECQUET	✓ Ready
Kit management		02/04/2025 11:35:19	SarcomaFusion		Julien Binet	✓ Ready
Product management		02/04/2025 11:34:12	LymphoSign		Julien Binet	✓ Ready
Configuration		01/04/2025 14:41:42	Sarcoma V2		Laurine BECQUET	✓ Ready
Help		17/03/2025 09:54:09	LymphoSign		Laurine BECQUET	✓ Ready
About		04/03/2025 14:38:50	LymphoSign		Laurine BECQUET	✓ Ready
Julien Binet jbinet@genexpath.com		17/02/2025 14:35:10	LymphoTranscript - RUO		Laurine BECQUET	✓ Ready
Genexpath Rouen genexpath		12/02/2025 10:18:44	Sarcoma		Laurine BECQUET	✓ Ready
Settings		05/02/2025 13:48:56	Sarcoma		Laurine BECQUET	✓ Ready
		Total Rows: 218				

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Back to Results

Ready

The following barcodes were not found

- GEP-BC-001
- GEP-BC-002
- GEP-BC-003
- GEP-BC-004

SEND REQUEST FOR REFUND OF BALANCE DEDUCTED

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RT-MIS		Results				
Home New analysis Results History Administration Analyses management History management Kit management Product management Configuration Help About		Date	Product	Internal reference	User	Status
Julien Binet jbinet@genexpath.com Genexpath Rouen Genexpath Settings		02/05/2025 11:29:21	LymphoSign		Laurine BECQUET	✓ Ready
		25/04/2025 15:07:40	CustomPanel		Laurine BECQUET	✓ Ready
		25/04/2025 15:02:02	LymphoTranscript - RUO		Laurine BECQUET	✓ Ready
		25/04/2025 14:32:27	SarcomaFusi		Julien Binet	✓ Ready
		25/04/2025 14:31:22	LymphoSign		Julien Binet	✓ Ready
		08/04/2025 10:55:16	Panel Huminflan		Laurine BECQUET	✓ Ready
		02/04/2025 11:35:19	SarcomaFusion		Julien Binet	✓ Ready
		02/04/2025 11:34:12	LymphoSign		Julien Binet	✓ Ready
		01/04/2025 14:41:42	Sarcoma V2		Laurine BECQUET	✓ Ready
		17/03/2025 09:54:09	LymphoSign		Laurine BECQUET	✓ Ready
		04/03/2025 14:38:50	LymphoSign		Laurine BECQUET	✓ Ready
		17/02/2025 14:35:10	LymphoTranscript - RUO		Laurine BECQUET	✓ Ready
		12/02/2025 10:18:44	Sarcoma		Laurine BECQUET	✓ Ready
		05/02/2025 13:48:56	Sarcoma		Laurine BECQUET	✓ Ready
						Total Rows: 218

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Status

✓ Ready

⚠ The following barcodes were not found

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