

Ne  *Biotech*

Mumps Virus IgG

NB-06-0389

Mumps Virus IgG

#Cat: NB-06-0389 Size: 1x96 Tests

1. INTRODUCTION

Mumps viruses are RNA viruses of the family Paramyxoviridae. The virions are spherical particles of 150-250 nm in diameter consisting of a ribonucleoprotein with helical symmetry and enveloped by matrix protein and a lipid bilayer which contains two spikeline structures: viral hemagglutinin (H) and viral neuraminidase (N). Mumps virus involves primarily the parotid and related salivary glands; however, infection can lead to CNS disease and accumulation of the virus in CSF. Mumps (Epidemic Parotitis) is an acute contagious viral disease mostly occurring in children. Nearly 50% of all infections are subclinical. The highest incidence of clinical manifestations is found in the age group of 4 to 15 years. Secondary infections are rare because of long-lasting immunity. 10 to 35 % of mumps cases develop orchitis which occurs nearly always after puberty. The process is mostly unilateral and the prognosis usually good. Mumps virus has been one of the most important causes of viral CNS disease (meningitis and encephalitis) in USA; vaccine administration has greatly reduced its incidence.

Species	Disease	Symptoms	Mechanism of Infection
Mumps virus	Mumps	Fever and unilateral or bilateral swelling of the parotid gland; the sublingual and submaxillary glands may also be involved Complications: Orchitis, Meningoencephalitis, Pancreatitis	Virus transmission occurs by droplet infection; invasion site is the upper respiratory tract. The incubation time period is 14 to 24 days

The presence of virus resp. infection may be identified by:

- Complement fixation, Hemadsorption, Hemagglutination-Inhibition
- Detection of antibody production by ELISA

Measurement of IgG/IgM antibodies in the serum allows for the determination of the immune status of the patients and recommendations about the necessity of immunization. As prophylaxis live attenuated virus is available and is frequently given in combination with others (measles-mumps-rubella). These vaccinations produce antibody formation in up to 90 % of vaccines.

2. INTENDED USE

The Neo Biotech Mumps Virus IgG-ELISA is intended for the qualitative determination of IgG class antibodies against Mumps virus in human serum or plasma (citrate).

3. PRINCIPLE OF THE ASSAY

The qualitative immunoenzymatic determination of IgG-class antibodies against Mumps Virus is based on the ELISA (Enzyme-linked Immunosorbent Assay) technique.

Microtiter strip wells are precoated with Mumps Virus antigens to bind corresponding antibodies of the specimen. After washing the wells to remove all unbound sample material horseradish peroxidase (HRP) labelled anti-human IgG conjugate is added. This conjugate binds to the captured Mumps Virus-specific antibodies. The immune complex formed by the bound conjugate is visualized by adding Tetramethylbenzidine (TMB) substrate which gives a blue reaction product. The intensity of this product is proportional to the amount of Mumps Virus-specific IgG antibodies in the specimen. Sulphuric acid is added to stop the reaction. This produces a yellow endpoint colour. Absorbance at 450 nm is read using an ELISA microwell plate reader.

4. MATERIALS

4.1. Reagents supplied

- **Mumps Virus Coated Wells (IgG):** 12 breakapart 8-well snap-off strips coated with Mumps Virus antigen; in resealable aluminium foil.
- **IgG Sample Diluent ***:** 1 bottle containing 100 ml of buffer for sample dilution; pH 7.2 \pm 0.2; coloured yellow; ready to use; white cap.
- **Stop Solution:** 1 bottle containing 15 ml sulphuric acid, 0.2 mol/l; ready to use; red cap.
- **Washing Solution (20x conc.) *:** 1 bottle containing 50 ml of a 20-fold concentrated buffer (pH 7.2 \pm 0.2) for washing the wells; white cap.
- **Mumps Virus anti-IgG Conjugate**:** 1 bottle containing 20 ml of peroxidase labelled rabbit antibody to human IgG; coloured blue, ready to use; black cap.
- **TMB Substrate Solution:** 1 bottle containing 15 ml 3,3',5,5'-tetramethylbenzidine (TMB); ready to use; yellow cap.
- **Mumps Virus IgG Positive Control***:** 1 bottle containing 2 ml; coloured yellow; ready to use; red cap.
- **Mumps Virus IgG Cut-off Control***:** 1 bottle containing 3 ml; coloured yellow; ready to use; green cap.
- **Mumps Virus IgG Negative Control***:** 1 bottle containing 2 ml; coloured yellow; ready to use; blue cap.

* Contains 0.1 % Bronidox L after dilution

** contains 0.2 % Bronidox L

*** contains 0.1 %

4.2. Materials supplied

- 1 Strip holder
- 1 Cover foil
- 1 Test protocol
- 1 distribution and identification plan

4.3. Materials and Equipment needed

- ELISA microwell plate reader, equipped for the measurement of absorbance at 450/620nm
- Incubator 37°C
- Manual or automatic equipment for rinsing wells
- Pipettes to deliver volumes between 10 and 1000 µl
- Vortex tube mixer
- Deionised or (freshly) distilled water
- Disposable tubes
- Timer

5. STABILITY AND STORAGE

The reagents are stable up to the expiry date stated on the label when stored at 2...8 °C.

6. REAGENT PREPARATION

It is very important to bring all reagents, samples and controls to room temperature (20...25°C) before starting the test run!

6.1. Coated snap-off Strips

The ready to use breakapart snap-off strips are coated with Mumps Virus antigen. Store at 2...8°C. *Immediately after removal of strips, the remaining strips should be resealed in the aluminium foil along with the desiccant supplied and stored at 2...8 °C; stability until expiry date.*

6.2. Mumps Virus anti-IgG Conjugate

The bottle contains 20 ml of a solution with anti-human-IgG horseradish peroxidase, buffer, stabilizers, preservatives and an inert blue dye. The solution is ready to use. Store at 2...8°C. *After first opening until expiry date when stored at 2...8°C.*

6.3. Controls

The bottles labelled with Positive, Cut-off and Negative Control contain a ready to use control solution. It contains 0.1% Kathon and has to be stored at 2...8°C. *After first opening until expiry date when stored at 2...8°C.*

6.4. IgG Sample Diluent

The bottle contains 100ml phosphate buffer, stabilizers, preservatives and an inert yellow dye. It is used for the dilution of the patient specimen. This ready to use solution has to be stored at 2...8°C. *After first opening until expiry date when stored at 2...8°C.*

6.5. Washing Solution (20xconc.)

The bottle contains 50 ml of a concentrated buffer, detergents and preservatives. Dilute washing solution 1+19; e.g. 10 ml washing solution + 190 ml fresh and germ-free redistilled water. The diluted buffer will keep for 5 days if stored at room temperature. Crystals in the solution disappear by warming up to 37 °C in a water bath. After first opening the concentrate is stable until the expiry date.

6.6. TMB Substrate Solution

The bottle contains 15 ml of a tetramethylbenzidine/hydrogen peroxide system. The reagent is ready to use and has to be stored at 2...8°C, away from the light. *The solution should be colourless or could have a slight blue tinge. If the substrate turns into blue, it may have become contaminated and should be thrown away. After first opening until expiry date when stored at 2...8°C.*

6.7. Stop Solution

The bottle contains 15 ml 0.2 M sulphuric acid solution (R 36/38, S 26). This ready to use solution has to be stored at 2...8°C. After first opening stability expiry date.

7. SPECIMEN COLLECTION AND PREPARATION

Use human serum or plasma (citrate) samples with this assay. If the assay is performed within 5 days after sample collection, the specimen should be kept at 2...8°C; otherwise, they should be aliquoted and stored deep-frozen (-20 to -70°C). If samples are stored frozen, mix thawed samples well before testing. *Avoid repeated freezing and thawing.*

Heat inactivation of samples is not recommended.

7.1. Sample Dilution

Before assaying, all samples should be diluted 1+100 with IgG Sample Diluent. Dispense 10µl sample and 1ml IgG Sample Diluent into tubes to obtain a 1+100 dilution and thoroughly mix with a Vortex.

8. ASSAY PROCEDURE

8.1. Test Preparation

Please read the test protocol carefully **before** performing the assay. Result reliability depends on strict adherence to the test protocol as described. The following test procedure is only validated for manual procedure. If performing the test on ELISA automatic systems we recommend to increase the washing steps from three to five and the volume of washing solution from 300µl to 350µl to avoid washing effects. Prior to commencing the assay, the distribution and identification plan for all specimens and controls should be carefully established on the result sheet supplied in the kit. Select the required number of microtiter strips or wells and insert them into the holder.

Please allocate at least:

1 well	(e.g. A1)	for the substrate blank,
1 well	(e.g. B1)	for the negative control,
2 wells	(e.g. C1+D1)	for the cut-off control and
1 well	(e.g. E1)	for the positive control.

It is recommended to determine controls and patient samples in duplicate, if necessary. Perform all assay steps in the order given and without any appreciable delays between the steps.

A clean, disposable tip should be used for dispensing each control and sample.

Adjust the incubator to $37^{\circ} \pm 1^{\circ}\text{C}$.

1. Dispense 100µl controls and diluted samples into their respective wells. Leave well A1 for substrate blank.
2. Cover wells with the foil supplied in the kit.
3. **Incubate for 1 hour \pm 5 min at $37\pm 1^{\circ}\text{C}$.**
4. When incubation has been completed, remove the foil, aspirate the content of the wells and wash each well three times with 300µl of Washing Solution. Avoid overflows from the reaction wells. The soak time between each wash cycle should be >5sec. At the end carefully remove remaining fluid by tapping strips on tissue paper prior to the next step!

Note: Washing is critical! Insufficient washing results in poor precision and falsely elevated absorbance values.

5. Dispense 100µl Mumps Virus anti-IgG Conjugate into all wells except for the blank well (e.g. A1). Cover with foil.
6. **Incubate for 30 min at room temperature.** *Do not expose to direct sunlight.*
7. Repeat step 4.
8. Dispense 100µl TMB Substrate Solution into all wells
9. **Incubate for exactly 15 min at room temperature in the dark.**
10. Dispense 100µl Stop Solution into all wells in the same order and at the same rate as for the TMB Substrate Solution.

Any blue colour developed during the incubation turns into yellow.

Note: Highly positive patient samples can cause dark precipitates of the chromogen!

These precipitates have an influence when reading the optical density. Predilution of the sample with physiological sodium chloride solution, for example 1+1, is recommended. Then dilute the sample 1+100 with dilution buffer and multiply the results in NTU by 2.

11. Measure the absorbance of the specimen at 450/620nm within 30 min after addition of

the Stop Solution.

8.2. Measurement

Adjust the ELISA Microwell Plate Reader **to zero** using the **substrate blank in well A1**.

If - due to technical reasons - the ELISA reader cannot be adjusted to zero using the substrate blank in well A1, subtract the absorbance value of well A1 from all other absorbance values measured in order to obtain reliable results!

Measure the absorbance of all wells at **450 nm** and record the absorbance values for each control and patient sample in the distribution and identification plan.

Dual wavelength reading using 620 nm as reference wavelength is recommended.

Where applicable calculate the **mean absorbance values** of all duplicates.

9. RESULTS

9.1. Run Validation Criteria

In order for an assay to be considered valid, the following criteria must be met:

- Substrate blank in A1 Absorbance value < **0.100**.
- Negative control in B1: Absorbance value < **0.200** and < **cut-off**
- Cut-off control in C1 and D1: Absorbance value **0.150 – 1.30**.
- Positive control in E1: Absorbance value > **cut-off**.

If these criteria are not met, the test is not valid and must be repeated.

9.2. Calculation of Results

The cut-off is the mean absorbance value of the Cut-off control determinations.

Example:

Absorbance value Cut-off control 0.39 + absorbance value Cut-off control 0.37 = 0.76 / 2 = 0.38

Cut-off = 0.38

9.3. Interpretation of Results

Samples are considered **POSITIVE** if the absorbance value is higher than 10% over the cut-off.

Samples with an absorbance value of 10% above or below the cut-off should not be considered as clearly positive or negative

→ **grey zone**

It is recommended to repeat the test again 2 - 4 weeks later with a fresh sample. If results in the second test are again in the grey zone

the sample has to be considered **NEGATIVE**.

Samples are considered **NEGATIVE** if the absorbance value is lower than 10% below the cut-off.

9.3.1. Results in Neo Biotech Units

$$\frac{\text{Patient (mean) absorbance value} \times 10}{\text{Cut-off}} = [\text{Neo Biotech-Units} = \text{NTU}]$$

Example:
$$\frac{1.786 \times 10}{0.38} = 47 \text{ NTU (Neo Biotech Units)}$$

Cut-off:	10	NTU
Grey zone:	9-11	NTU
Negative:	<9	NTU
Positive:	>11	NTU

10. SPECIFIC PERFORMANCE CHARACTERISTICS

10.1. Precision

<u>Interassay</u>	<u>n</u>	<u>Mean</u>	<u>Cv (%)</u>
Pos. Serum	12	28.1	2.9
<u>Intraassay</u>	<u>n</u>	<u>Mean</u>	<u>Cv (%)</u>
Pos. Serum	7	1.10	5.4

10.2. Diagnostic Specificity

The diagnostic specificity is defined as the probability of the assay of scoring negative in the absence of the specific analyte. It is 87.5 %.

10.3. Diagnostic Sensitivity

The diagnostic sensitivity is defined as the probability of the assay of scoring positive in the presence of the specific analyte. It is >95 %.

10.4. Interferences

Interferences with hemolytic, lipemic or icteric sera are not observed up to a concentration of 10 mg/ml hemoglobin, 5 mg/ml triglycerides and 0.2 mg/ml bilirubin.

Note: The results refer to the groups of samples investigated; these are not guaranteed specifications.

11. LIMITATIONS OF THE PROCEDURE

Bacterial contamination or repeated freeze-thaw cycles of the specimen may affect the absorbance values. Diagnosis of an infectious disease should not be established on the basis of a single test result. A precise diagnosis should take into consideration clinical history, symptomatology as well as serological data.

In immunocompromised patients and newborns serological data only have restricted value.

12. PRECAUTIONS AND WARNINGS

- In compliance with article 1 paragraph 2b European directive 98/79/EC the use of the

in vitro diagnostic medical devices is intended by the manufacturer to secure suitability, performances and safety of the product. Therefore, the test procedure, the information, the precautions and warnings in the instructions for use have to be strictly followed. The use of the test kits with analyzers and similar equipment has to be validated. Any change in design, composition and test procedure as well as for any use in combination with other products not approved by the manufacturer is not authorized; the user himself is responsible for such changes. The manufacturer is not liable for false results and incidents for these reasons. The manufacturer is not liable for any results by visual analysis of the patient samples.

- Only for in-vitro diagnostic use.
- All components of human origin used for the production of these reagents have been tested for anti-HIV antibodies, anti-HCV antibodies and HBsAg and have been found to be non-reactive. Nevertheless, all materials should still be regarded and handled as potentially infectious.
- Do not interchange reagents or strips of different production lots.
- No reagents of other manufacturers should be used along with reagents of this test kit.
- Do not use reagents after expiry date stated on the label.
- Use only clean pipette tips, dispensers, and lab ware.
- Do not interchange screw caps of reagent vials to avoid cross-contamination.
- Close reagent vials tightly immediately after use to avoid evaporation and microbial contamination.
- After first opening and subsequent storage check conjugate and control vials for microbial contamination prior to further use
- To avoid cross-contamination and falsely elevated results pipette patient samples and dispense conjugate without splashing accurately to the bottom of wells.
- The Neo Biotech ELISA is only designed for qualified personnel who are familiar with good laboratory practice.

WARNING: In the used concentration Bronidox L has hardly any toxicological risk upon contact with skin and mucous membranes!

WARNING: Sulphuric acid irritates eyes and skin. Keep out of the reach of children. Upon contact with the eyes, rinse thoroughly with water and consult a doctor!

12.1. Disposal Considerations

Residues of chemicals and preparations are generally considered as hazardous waste. The disposal of this kind of waste is regulated through national and regional laws and regulations. Contact your local authorities or waste management companies which will give advice on how to dispose hazardous waste.

13. ORDERING INFORMATION

Prod. No.: NB-06-0389 Mumps Virus IgG-ELISA (96 Determinations)





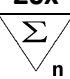
BIBLIOGRAPHY

Berbers, G.A.M et al.: Blocking ELISA for detection of mumps virus antibodies in human sera. J. Virol. Methods 42 (1993) 155-168
Berger, R.,M. Just: Comparison of five different tests for mumps antibodies. Infection 8 (1980) 180-183

Forsey, T.,M.L. Bentley, P.D. Minor: Mumps vaccines and meningitis. Lancet 340 (1992) 980

Friedman, M., I. Hadari et al.: Virus-specific secretory IgA antibodies as a means of rapid diagnosis of measles and mumps infection. Isr. J. Med. Sci. 19 (1993) 881-883

Gerike, E., A. Tischer: Erfahrung mit der Schutzimpfung gegen Masern, Mumps und Röteln in den neuen Bundesländern. Gesundheitswesen 55 (1993) 38-29

Symbols Key	
	Manufactured by
IVD	In Vitro Diagnostic Medical Device
LOT	Lot Number
	Expiration Date
	Storage Temperature
CE	CE Mark
[REF]	Catalogue Number
	Consult Instructions for Use
MTP	Microplate
CONJ	Conjugate
CONTROL -	Control serum, negative
CONTROL +	Control serum, positive
CUT OFF	Cut off control serum
DIL G	Sample diluent buffer IgG
SOLN STOP	Stop solution
SUB TMB	TMB Substrate solution
WASHBUF 20x	Washing solution 20x concentrated
	Contains sufficient for "n" tests

SCHEME OF THE ASSAY

Mumps Virus IgG-ELISA

Test Preparation

Prepare reagents and samples as described.
 Establish the distribution and identification plan for all specimens and controls on the result sheet supplied in the kit.
 Select the required number of microtiter strips or wells and insert them into the holder.

Assay Procedure

	Substrate blank (e.g. A1)	Negative control	Positive control	Cut-off control	Sample (diluted 1+100)
Negative control	-	100µl	-	-	-
Positive control	-	-	100µl	-	-
Cut-off control	-	-	-	100µl	-
Sample (diluted 1+100)	-	-	-	-	100µl
Cover wells with foil supplied in the kit Incubate for 1 h at 37°C Wash each well three times with 300µl of washing solution					
Conjugate	-	100µl	100µl	100µl	100µl
Cover wells with foil supplied in the kit Incubate for 30 min at room temperature Wash each well three times with 300µl of washing solution					
TMB Substrate	100µl	100µl	100µl	100µl	100µl
Incubate for exactly 15 min at room temperature in the dark					
Stop Solution	100µl	100µl	100µl	100µl	100µl
Photometric measurement at 450 nm (reference wavelength: 620 nm)					

engl10082009-JH