

Technical Data Sheet

PE Mouse anti-eNOS

Product Information


Material Number:	560103
Alternate Name:	NOS type III, NOS3, EC-NOS, NOS III
Size:	50 tests
Vol. per Test:	20 µl
Clone:	33/eNOS
Immunogen:	Human eNOS aa. 1025-1203
Isotype:	Mouse IgG1, κ
Reactivity:	Confirmed: Human Reported: Mouse
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

Nitric oxide synthase (NOS), a cell-type specific enzyme, catalyzes the synthesis of nitric oxide (NO). NO is a short-lived radical that transmits signals involved in vasorelaxation, neurotransmission, and cytotoxicity. In neurons and endothelial cells, constitutive NOS (cNOS) is activated by agonists that increase intracellular Ca²⁺ levels and enhance calmodulin binding. Neuronal NOS (nNOS) and endothelial NOS (eNOS) have recognition sites for NADPH, FAD, FMN, and calmodulin and both are regulated in a similar manner. The human forms exhibit 52% amino acid identity. However, they are distinct gene products of about 155 kDa (nNOS) and 140 kDa (eNOS). The eNOS gene was cloned from human vascular endothelium as well as from bovine aortic endothelial cells (BAEC). eNOS protein has a unique N-myristylation consensus sequence that may explain its membrane localization.

The 33/eNOS monoclonal antibody recognizes eNOS, regardless of phosphorylation status.

The specificity of this antibody conjugate for flow cytometric analysis was validated by confirming that RNA-mediated interference (RNAi) of the specific protein reduced the staining of the cells (see figure). Furthermore, the capacity of the RNAi to down-regulate the expression of the relevant protein was confirmed by western blot analysis.

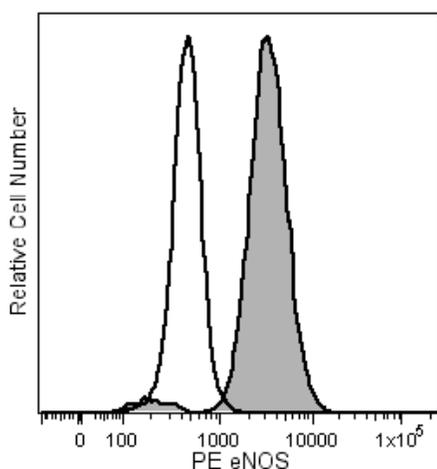

Analysis of eNOS in human endothelial cells. EA-hy 926 cells (Edgell, McDonald, Graham, 1983) were either transfected with eNOS RNAi (open histogram) or untreated (shaded histogram). The cells were fixed (BD Cytotfix™ Fixation buffer, Cat. No. 554655) for 10 minutes at 37°C, then permeabilized (BD™ Phosflow Perm Buffer III, Cat. No. 558050) on ice for 30 minutes, and then stained with PE Mouse anti-eNOS. Down-regulation of eNOS expression is evident in the RNAi-transfected cells. Flow cytometry was performed on a BD™ LSR II flow cytometry system.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with R-PE under optimum conditions, and unconjugated antibody and free PE were removed.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.268.5430	32.53.720.550	0120.8555.90	65.6861.0633	0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



Application Notes

Application

Intracellular staining (flow cytometry)

Routinely Tested

Recommended Assay Procedure:

Either BD Cytotfix™ fixation buffer or BD™ Phosflow Fix Buffer I may be used for cell fixation.

Suggested Companion Products

<u>Catalog Number</u>	<u>Name</u>	<u>Size</u>	<u>Clone</u>
558050	Perm Buffer III	125 ml	(none)
554655	Fixation Buffer	100 ml	(none)
557870	Fix Buffer I	250 ml	(none)
559320	PE Mouse IgG1, κ Isotype Control	100 tests	MOPC-21

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
3. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
4. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at wwwbdbiosciences.com/colors.
5. Please refer to wwwbdbiosciences.com/pharming/protocols for technical protocols.

References

- Chen PF, Tsai AL, Wu KK. Cysteine 184 of endothelial nitric oxide synthase is involved in heme coordination and catalytic activity. *J Biol Chem.* 1994; 269(40):25062-25066. (Clone-specific: Western blot)
- Dinerman JL, Dawson TM, Schell MJ, Snowman A, Snyder SH. Endothelial nitric oxide synthase localized to hippocampal pyramidal cells: implications for synaptic plasticity. *Proc Natl Acad Sci U S A.* 1994; 91(10):4214-4218. (Biology)
- Edgell C-JS, McDonald CC, Graham JB. Permanent cell line expressing human factor VIII-related antigen established by hybridization. *Proc Natl Acad Sci U S A.* 1983; 80:3734-3737. (Methodology: Controls)
- Nathan C, Xie QW. Regulation of biosynthesis of nitric oxide. *J Biol Chem.* 1994; 269(19):13725-13728. (Biology)
- Shen YH, Zhang L, Utama B et al. Human cytomegalovirus inhibits Akt-mediated eNOS activation through upregulating PTEN (phosphatase and tensin homolog deleted on chromosome 10). *Cardiovasc Res.* 2006; 69(2):502-511. (Biology)
- Varghese P, Harrison RW, Lofthouse RA, Georgakopoulos D, Berkowitz DE, Hare JM. β 3-adrenoceptor deficiency blocks nitric oxide-dependent inhibition of myocardial contractility. *J Clin Invest.* 2000; 106(5):697-703. (Clone-specific: Western blot)